BRAUN

IRT 3020 CO

Type 6026


ThermoScan

English	3
Français	9
Español	16

Internet:
www.kaz.com/braun

Distributed by:
Kaz, Incorporated
Hudson,
New York 12534
USA

6-026-360/00/VIII-07/G3
USA/CDN/MEX

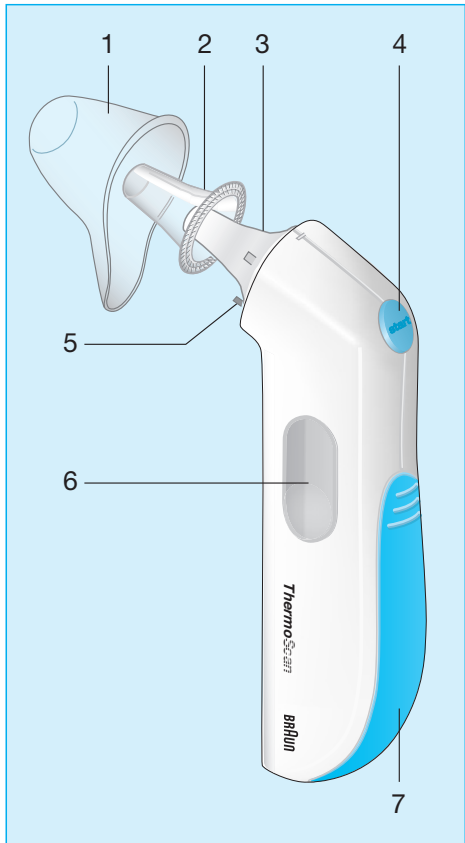English

The Braun ThermoScan thermometer has been carefully developed for accurate, safe and fast temperature measurements in the ear. The shape of the thermometer probe prevents it from being inserted too far into the ear canal and hurting the eardrum.

However, as with any thermometer, proper technique is very important to obtaining accurate temperatures. Therefore, read the use instructions carefully and thoroughly.

Important

- The operating ambient temperature range for this thermometer is 50–104 °F (10–40 °C).
- Do not expose the thermometer to temperature extremes (below –4 °F / –20 °C or over 122 °F / 50 °C) nor excessive humidity (> 95 % RH).
- This thermometer must only be used with genuine Braun ThermoScan Lens Filters. Other lens filters can lead to inaccuracy.
- To avoid inaccurate measurements always use this thermometer with a new, clean lens filter attached.
- If the thermometer is accidentally used without a lens filter attached, clean the lens (see «Care and cleaning»).
- Keep lens filters out of reach of children.
- This thermometer is intended for household use only.
- Use of this thermometer is not intended as a substitute for consultation with your physician.

Body temperature

Normal body temperature is a range. The following table shows that this normal range also varies by site. Therefore, measurements from different sites should not be directly compared.

Normal ranges by site:

Axillary ¹ :	94.5 – 99.1 °F	34.7 – 37.3 °C
Oral ¹ :	95.9 – 99.5 °F	35.5 – 37.5 °C
Rectal ¹ :	97.9 – 100.4 °F	36.6 – 38.0 °C
ThermoScan ² :	96.4 – 100.4 °F	35.8 – 38.0 °C

- A person's normal temperature range tends to decrease with age. The following table shows normal ThermoScan ranges by age.

Normal ThermoScan ranges by age²:

0 – 2 years	97.5 – 100.4 °F	36.4 – 38.0 °C
3 – 10 years	97.0 – 100.0 °F	36.1 – 37.8 °C
11 – 65 years	96.6 – 99.7 °F	35.9 – 37.6 °C
> 65 years	96.4 – 99.5 °F	35.8 – 37.5 °C

The range of normal also varies from person to person and fluctuates throughout the day. It is therefore important to determine normal temperature ranges. This is easily done using Braun ThermoScan. Practice taking temperatures on yourself and healthy family members to determine each person's normal temperature range.


Note: When consulting your physician, communicate that the ThermoScan temperature is a temperature measured in the ear and if possible, note the individual's normal ThermoScan temperature range as additional reference.

1. Chamberlain, J.M., Terndrup, T.E., New Light on Ear Thermometer Readings, Contemporary Pediatrics, March 1994.
2. Chamberlain, J.M., et al., Determination of Normal Ear Temperature with an Infrared Emission Detection Thermometer, Annals of Emergency Medicine, January 1995.

Product description

- 1 Protective cap
- 2 Lens filter
- 3 Probe
- 4 «start» button
- 5 Lens filter detector
- 6 Display
- 7 Battery door

How to use your Braun ThermoScan

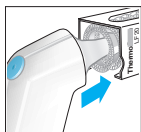1. To achieve accurate measurements, make sure a new, clean lens filter (2) is in place before each measurement.

2. To turn the thermometer on, push the «start» button (4).



During an internal self-check, the display (6) shows all segments. Then the last temperature taken will be displayed together with «MEM».


As soon as this ready symbol appears, you can start the measuring process.

3. Perform an ear tug to straighten the ear canal. This gives the thermometer a clear view of the eardrum.

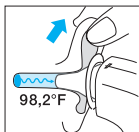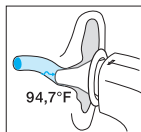
An ear tug is best performed by using your free hand to grasp the outer edge of the top half of the ear. To take your own temperature, wrap your free hand around the back of your head and tug your ear from behind.

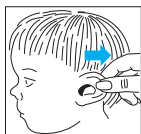Children under 1 year:
Pull the ear straight back.

Children aged 1 year to adult:
Pull the ear up and back.

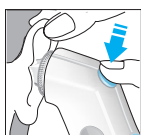4. While tugging the ear, fit the probe (3) snugly into the ear canal, then push the «start» button. A short beep will signal the start of the measuring process.

After one second a longer beep will signal the end of the measuring process. The result is shown on the display (6).

5. For the next measurement, remove the used lens filter and put on a new, clean lens filter. Clear the display by pushing the «start» button once. Fit the probe snugly into the ear canal, then push the «start» button.


The Braun ThermoScan ear thermometer turns off automatically after 60 seconds of inactivity.

Temperature taking hints

- The right ear measurement may differ from the measurement taken in the left ear. Therefore, always take the temperature in the same ear.
- The ear must be free from obstructions or excess earwax build-up to take an accurate reading.
- External factors may influence ear temperatures, including when an individual has:
 - been lying on one ear or the other
 - had their ears covered
 - been exposed to very hot or very cold temperatures, or
 - been recently swimming or bathing.
 In these cases, remove the individual from the situation and wait 20 minutes prior to taking a temperature.
- Use the untreated ear if ear drops or other ear medications have been placed in the ear canal.

Care and cleaning

The probe tip is the most delicate part of the thermometer. It has to be clean and intact to ensure accurate readings.

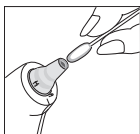If the thermometer is ever accidentally used without a lens filter, clean the probe tip as follows:
Very gently wipe the surface with a cotton swab or soft cloth moistened with alcohol.

After the probe tip has completely dried, you can put a new lens filter on and take a temperature measurement. If the probe tip is damaged, contact Kaz.

Use a soft, dry cloth to clean the thermometer display and exterior. Do not use abrasive cleaners. Never submerge this thermometer in water or any other liquid.

Store thermometer and lens filters in a dry location free from dust and contamination and away from direct sunlight.

Additional lens filters (LF 40) are available at most stores carrying Braun ThermoScan.

Trouble-shooting

Situation	Solution
 No lens filter is attached.	Attach new, clean lens filter.
 Ambient temperature is not within the allowed operating range (50–104 °F or 10–40 °C).	Allow the thermometer to remain for 30 minutes in a room where the temperature is between 50 and 104 °F or 10 and 40 °C .

 HI = too high  LO = too low Temperature taken is not within typical human temperature range (93.2–108 °F or 34 – 42.2 °C).	Make sure the probe tip and lens are clean and a new, clean lens filter is attached. Make sure the thermometer is properly inserted into the ear. Then, take a new temperature measurement.
 System error If error persists, If error still persists,	Wait 1 minute until the thermometer turns off automatically, then turn on again. ... reset the thermometer by removing the battery and putting it back in. ... call Kaz Customer Service at 1-800-327-7226
 Battery is low, but thermometer will still operate correctly.	Insert new battery.
 Battery is too low to take correct temperature measurement.	Insert new battery.
Do you have any further questions?	Call Kaz Customer Service.

Changing the temperature scale


Your Braun ThermoScan is shipped with the Fahrenheit (°F) temperature scale activated. If you wish to switch to Celsius (°C), proceed as follows:

Make sure the thermometer is turned off. Press and hold down the «start» button. After about 8 seconds the display will show this sequence: «°F» / «°C» /...


Release the «start» button when «°C» is shown. There will be a short beep to confirm the new setting, then the thermometer is turned off automatically.

Replacing the battery

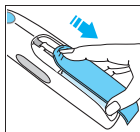The thermometer is supplied with one 1.5 V type AAA (LR3). For best performance, we recommend Duracell® alkaline batteries. Insert new battery when the battery symbol appears on the display.

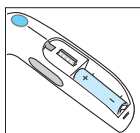Open the battery compartment. Remove the battery and replace with new battery, making sure the poles are in the right direction. Slide battery door back until it snaps in place.

To protect the environment, dispose of empty batteries at appropriate collection sites according to national or local regulations.

Calibration

The thermometer is initially calibrated at the time of manufacture. If this thermometer is used according to the use instructions, periodic re-adjustment is not required. If at any time you question the accuracy of temperature measurements, please contact Kaz Customer Service at 1-800-327-7226.

Manufacturing date is given by the LOT number located in the battery compartment. The first number after LOT represents the last digit of the year of manufacture.

The last two digits give the week of manufacture. An example: LOT 503 – this product has been manufactured in week 03 of the year 2005.

Product specifications

Displayed temperature range: 93.2 °F – 108 °F
(34 °C – 42.2 °C)

Operating ambient temperature range: 50 °F – 104 °F
(10 °C – 40 °C)

Display resolution: 0.1 °F or °C

Accuracy for patient temperature range

96.8 °F – 102.2 °F

(36 °C – 39 °C):

outside this range:

Maximum Laboratory Error

± 0.4 °F (± 0.2 °C)

± 0.5 °F (± 0.3 °C)

Long term storage ranges

Temperature: –4 °F to 122 °F
(–20 °C to 50 °C)

Humidity: ≤ 95 % non-condensing

Battery life: 2 years / 1000 measurements

If device is not used within specified temperature and humidity ranges the technical accuracy can not be ensured.

This infrared thermometer meets requirements established in ASTM Standard E 1965-98 (for the thermometer system [thermometer with lens filter]). Full responsibility for the conformance of the product to the standard is assumed by Kaz, Inc., Hudson, NY 12534 USA.

ASTM laboratory accuracy requirements for the thermometer only (not including a lens filter) in the display range of 96.8 °F to 102.2 °F (36 °C to 39 °C) for infrared thermometers is ± 0.4 °F (± 0.2 °C), whereas for mercury-in-glass and electronic thermometers, the requirement per ASTM Standards E 667-86 and E 1112-86 is ± 0.2 °F (± 0.1 °C).

This appliance conforms to the following standards:

DIN EN 60601-1: 3/96 «Medical electrical equipment» – Part 1: General requirements for safety

DIN EN 60601-1-2/2001 «Medical electrical equipment» – Part 1-2: General requirements for safety – collateral standard: Electromagnetic compatibility – requirements and tests

DIN EN 12470-5: 2003 «Clinical thermometers» – Part 5: Performance of infrared ear thermometers (with maximum device)

TYPE BF
EQUIPMENT

Attention, consult
ACCOMPANYING
DOCUMENTS

Internally Powered Equipment
Continuous Operation
Not Protected against Ingress of Water
U.S. Patent No. 5,088,834 Other Patents Pending

MEDICAL ELECTRICAL EQUIPMENT needs special precautions regarding EMC. For detailed description of EMC requirements please contact the Kaz customer service at 1-800-327-7226.

Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT.

Limited three year warranty

Kaz warrants this product against any defects that are due to faulty material or workmanship for a period of three years from the original date of consumer purchase or receipt as a gift. This warranty applies when used for normal household use in accordance with the Use and Care Manual and excludes the battery and damage to the product resulting from accident or misuse. This product is not warranted when used in a professional environment.

In no event shall Kaz be liable for any special, incidental, indirect or consequential damages in connection with the purchase or use of this product or costs over the original cost of the product.

If the product should not perform to specifications within the warranty period, contact Kaz Customer Support at 1 (800) 327-7226 for repair or replacement at no charge. Kaz reserves the right to replace a non conforming product with the most comparable product currently available. This product has a five year minimum expected lifetime.

While this warranty gives you specific legal rights, you may also have other rights which vary from state to state/province to province.

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